


FC161 Hand-held Fundus Camera




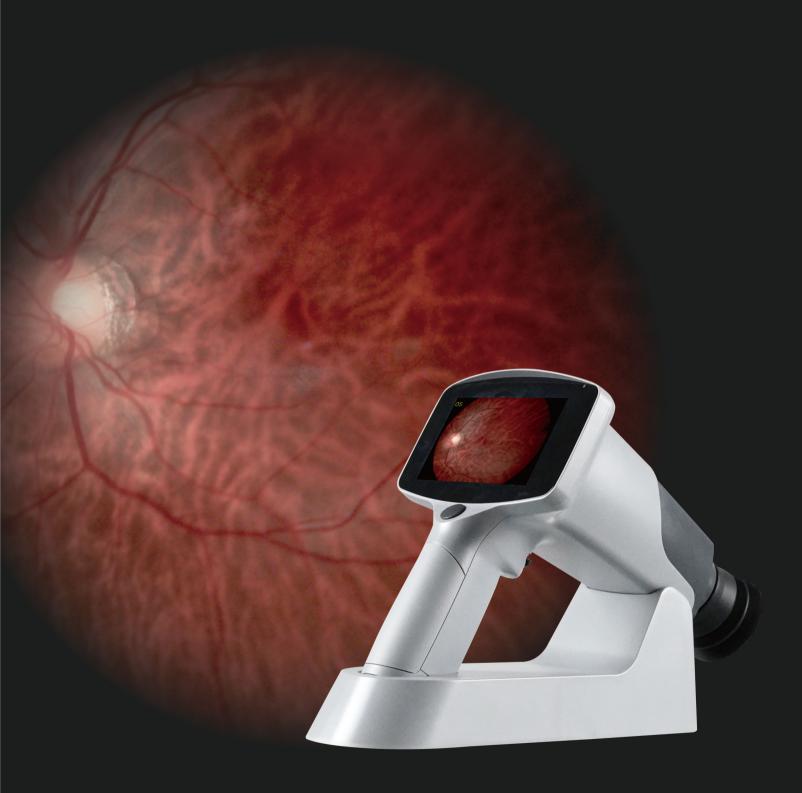


Auto Split Focusing Technology Millisecond Focusing Speed

FC161 adopts Auto Split Focusing technology, focusing speed reaches milliseconds.

It avoids the out-of-focus of the fundus image caused by the patient's movement during shooting, effectively improve the image quality and reduce the requirements for cooperation between operators and patients.

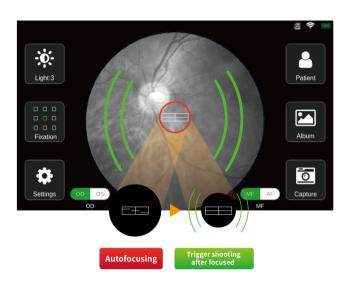

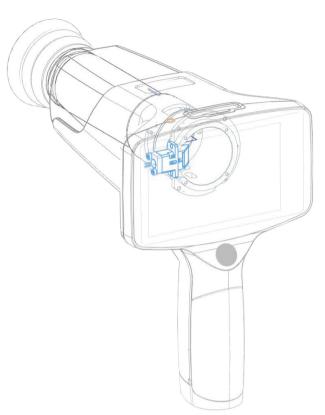


Advantages

High Resolution

FC161 provides high resolution fundus images,image resolution is up to 12 Mega Pixels.

It makes the diagnosis of early lesions more accurate and also essential for diagnosing retina disease, such as diabetic retinopathy,glaucoma and macular degeneration.



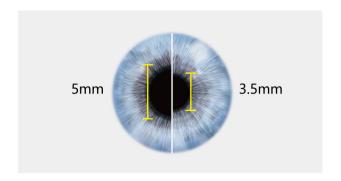


Non-Mydriatic

The minimum pupil diameter of FC161 is ϕ 3.5mm. 20 minutes of mydriatic waiting time is saved, and the risk of mydriasis to the patient is avoided.


9 Fixation Targets

With 9 fixation targets, FC161 covers 85° field of view of the fundus and supports the examination of early peripheral lesions of the fundus.

With $45\ ^\circ$ field of view , a single central fundus image can be taken to meet the needs of fundus disease screening.

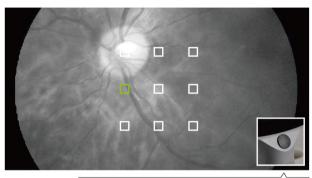

Easily select different fixation target by rotating direction key.

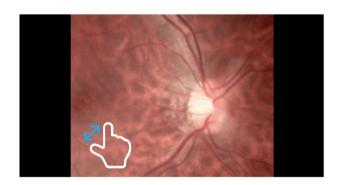

Smart Touch Operation

4.3 " full-touch LED screen brings a smarter operating experience. Doctors can zoom the image at any time to check the details of the fundus, and can also easily slide the screen to edit and view the case.

Real-time Uploading

FC161 supports real-time synchronization between fundus images and Mediview Patient Management software through WIFI connection.



Specification

FC161
Non-mydriatic
45degrees
12M Pixels
9×LED,Internal
3.5mm
White LED & IR LED
Manual/Automatic
-20D ~ +20D
4.3" Touch Screen
JPEG
Color
Micro SD card. Maximum 32G
WIFI/USB
Rechargeable battery 3.7V/3400mAh x 2
3 hours of continuous work
100-240V~, 0.5A, 50/60Hz
5V DC,2A
28cm×13cm×15cm
800g
Chinese, English

Standard Accessories





Hand-held Fundus Camera

USB

Cable



Recharging Station



Protection Cover



Power Adaptor







SD Card



Mediview Patient Management Software

Optional Accessories



Slit Lamp Attachment



(€ ₽ NMPA

Shanghai MediWorks Precision Instruments Co.,Ltd.

Add:Building 7, Ming Pu Plaza, No. 3279, San Lu Rd, Min Hang District, Shanghai, 201100, China Tel: +86-21-54260421 54260423 Fax:+86-21-54260425 Email: Marketing@mediworks.biz International@mediworks.biz

Follow Us on







www.mediworks.biz